




®

TERUMO SYRINGE

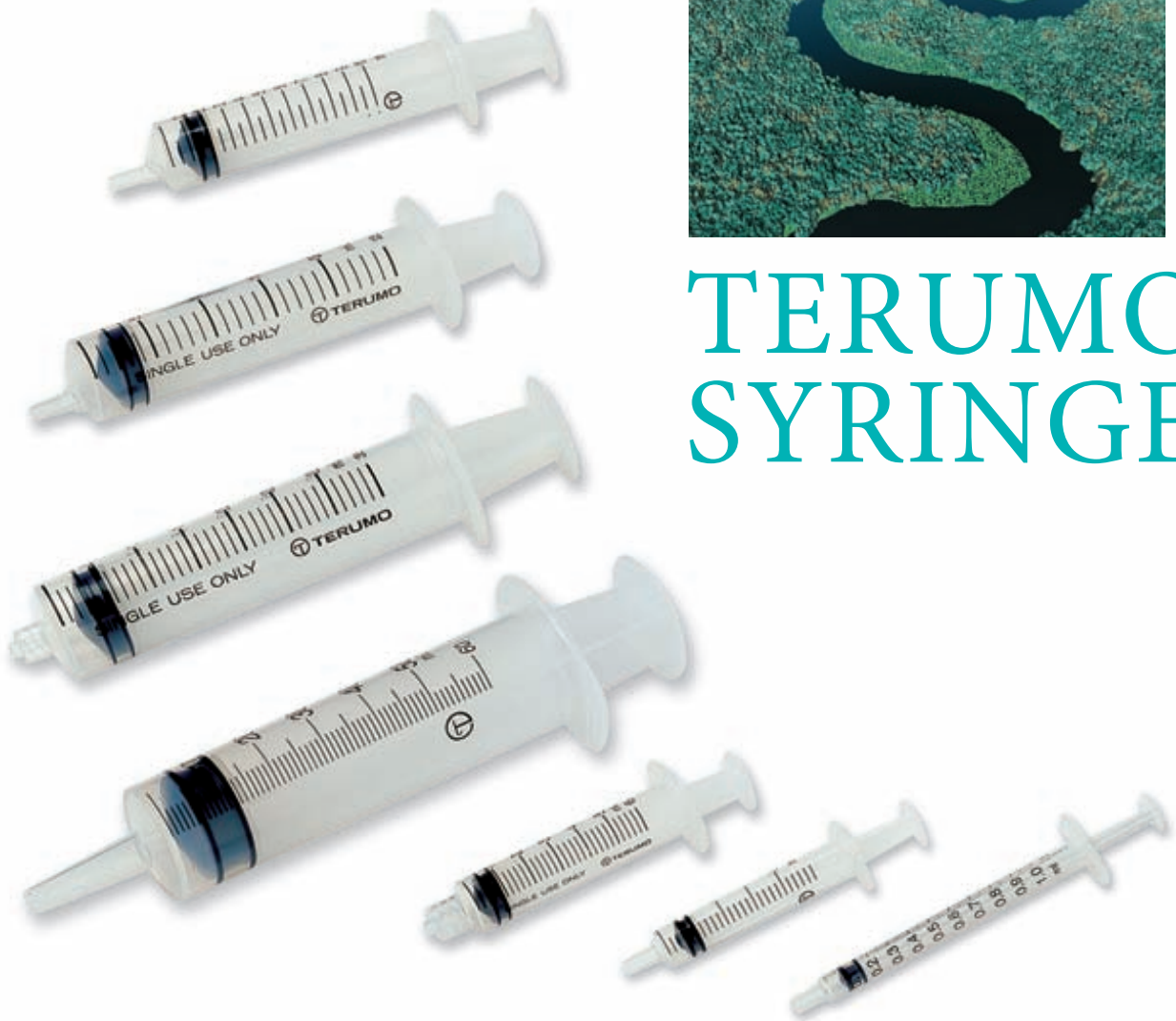

Terumo defines a new era of global quality

- Double-contact gasket for an airtight seal
- Advanced siliconisation for optimal plunger movement
- Firm and definite end-stop
- Clear scale markings to enable precise dosage

- Highly transparent barrel for visualisation
- Ergonomic finger grips
- Latex and PVC-free



A smooth and perfectly controlled injection is possible due to the double-contact gasket to ensure an airtight seal. It makes a smooth and precise plunger action possible.

The ergonomically designed finger grips and push-button allow a perfectly controlled and stable injection.



The highly transparent barrel allows for easy recognition of air bubbles and turbid medication. The bold scale markings on the barrel assure the precise measurement of the medication dose.

The firm and definite end stop prevents accidental separation or spillages, and avoids exposure to hazardous or infectious fluids.

We equally offer an OPAQUE syringe within our range that can be used to administer photosensitive drugs. The syringe is latex free, opaque, 50 ml, with a luer lock system and can be used on syringe pumps.

The latex free gasket designed for safety uses a highly inert thermoplastic elastomer material. It reduces the risk of leaking or medication interaction. This unique material equally allows the plunger to move more smoothly.

Our sterile, individual blister packs keep every syringe in a sterile condition until it will be used. It also provides tamper-proof evidence and is latex free.



Syringes without needle

Product Code	Description		Packaging	
	Size (ml)	Type	Unit box	Shipping carton
BS-01T	1	Luer Slip tip	100	1800
SS-02S	2	Luer Slip tip	100	1200
SS*02LE1*	2.5	Luer Lock tip	100	2400
SS+03L*	3	Luer Lock tip	100	1800
SS-05S	5	Luer Slip tip	100	1000
SS*05LE1*	5	Luer Lock tip	100	1600
SS+05L1*	5	Luer Lock tip	100	1200
SS+10ES	10	Luer Eccentric tip	100	500
SS+10L*	10	Luer Lock tip	100	800
SS+20ES*	20	Luer Eccentric tip	50	600
SS+20L*	20	Luer Lock tip	50	600
SS*30ESE1*	30	Luer Eccentric tip	50	400
SS*30LE1*	30	Luer Lock tip	50	400
BS-50ES	50	Luer Eccentric tip	25	100
BS-50LG	50	Luer Lock tip	25	100
BS-50C	50	Luer Catheter tip	25	100
BS-50LB	50	Luer Opaque Lock tip	25	100
SS*50LE*	50	Luer Lock tip	20	240
SS*50ESE*	50	Luer Slip tip	20	240

Syringes with needle

Product Code	Description			Packaging	
	Size (ml)	Diameter cannula	Length cannula	Unit box	Shipping carton
BS-01T2516	1	25G (0.5 mm)	5/8" (16 mm)	100	1800
BS-01T2613	1	26G (0.45 mm)	1/2" (12 mm)	100	1800
SS-02S2116	2	21G (0.8 mm)	5/8" (16 mm)	100	600
SS-02S2125	2	21G (0.8 mm)	1" (25 mm)	100	600
SS-02S2138	2	21G (0.8 mm)	1 1/2" (40 mm)	100	600
SS-02S2238	2	22G (0.7 mm)	1 1/2" (40 mm)	100	600
SS-02S2325	2	23G (0.6 mm)	1" (25 mm)	100	600
SS-02S2332	2	23G (0.6 mm)	1 1/4" (30 mm)	100	600
SS-02S2516	2	25G (0.5 mm)	5/8" (16 mm)	100	600
SS-05S2138	5	21G (0.8 mm)	1 1/2" (40 mm)	100	800
SS-05S2238	5	22G (0.7 mm)	1 1/2" (40 mm)	100	800
SS-05S2332	5	23G (0.6 mm)	1 1/4" (30 mm)	100	800
SS-10S2038	10	20G (0.9 mm)	1 1/2" (40 mm)	100	600
SS-10S2138	10	21G (0.8 mm)	1 1/2" (40 mm)	100	600
SS-10S2238	10	22G (0.7 mm)	1 1/2" (40 mm)	100	600

Sterilized by ethylene oxide
* Sterilized by e-beam



Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, **Terumo Europe** provides you with the clinical systems, means and solutions to keep life flowing.



TERUMO CORPORATION
44-1, 2-chome Hatagaya Shibuya-ku Tokyo 151-0072 JAPAN
Tel. +81/3 3374.8111 Fax +81/3 3374.8196

TERUMO EUROPE N.V.
Researchpark Haasrode 1520 Interleuvenlaan 40 B-3001 Leuven BELGIUM
Tel. +32/16.38.12.11 Fax +32/16.40.02.49

TERUMO EUROPE N.V. BENELUX SALES DIVISION
Researchpark Haasrode 1515 Interleuvenlaan 40 B-3001 Leuven BELGIUM
Tel. +32/16.39.25.80 Fax +32/16.39.25.99 THE NETHERLANDS
Tel. 0800/022.03.96 Fax 0800/022.04.14

LABORATOIRES TERUMO FRANCE S.A.
1, Parc Ariane Bâtiment Uranus Rue Hélène Boucher 78284 Guyancourt Cedex FRANCE
Tel. +33/1 30.96.13.00 Fax +33/1 30.43.60.85

TERUMO DEUTSCHLAND GmbH
Hauptstrasse 87 65760 Eschborn GERMANY
Tel. +49/61.96.80.230 Fax +49/61.96.80.23.200

TERUMO DEUTSCHLAND GmbH
Zweigniederlassung Wettingen Bahnhofstrasse 88 CH-5430 Wettingen SWITZERLAND
Tel. +41/56.419.10.10 Fax +41/56.419.10.11

TERUMO EUROPE ESPAÑA SL
Avda. Juan Carlos I, N° 13 - 7ª Planta Edificio Torre La Garena 28806 Alcalá de Henares (Madrid) SPAIN
Tel. general: +34/902.10.12.98
Tel. cust. serv.: +34/902.20.31.40
Fax general: +34/902.10.13.58
Fax cust. serv.: +34/902.10.17.57

TERUMO ITALIA SRL
Via Simone Martini 143/145 I-00142 Roma ITALY
Tel. +39/06.519.61.420 Fax +39/06.503.04.07

TERUMO SWEDEN AB
Nya Varvet, Byggnad 90 Sven Källfelts Gata 18 SE-426 71 Västra Frölunda SWEDEN
Tel. +46/31.748.58.80 Fax +46/31.748.58.90

TERUMO DENMARK Filial af TERUMO SWEDEN AB
Tamesis The Causeway Egham Surrey TW20 9AW UNITED KINGDOM
Tel. +44/1784.476.200 Fax +44/1784.476.201

TERUMO UK Ltd.
Tamesis The Causeway Egham Surrey TW20 9AW UNITED KINGDOM
Tel. +44/1784.476.200 Fax +44/1784.476.201

TERUMO EUROPE N.V.
Moscow Representative office 4, Ivana Franko St 121108 Moscow RUSSIA
Tel. +7/495.234.3303 Fax. +7/495.234.3303

©: Registered Trademark Published by Terumo Europe N.V.

Terumo Vietnam Medical Equipment Co., Ltd. Head office
14th & 16th Floor, Geleximco Building, 36 Hoang Cau Street, O Cho Dua Ward, Dong Da District Hanoi, Vietnam
Tel: +84 24 3936 1643 Fax: +84 24 3936 1641

Terumo Vietnam Medical Equipment Co., Ltd. HCMC Branch office
Room 505, Hoang Anh Safomec Building, 7/1 Thanh Thai, Ward 14, District 10 HCM City, Vietnam
Tel: +84 28 3866 9263

www.terumo.com
www.terumo-europe.com

TVME-2310-004